


Small Animal Anesthesia Machines ▶

R540 Enhanced Anesthesia Machine

Regular assembly (R540-A), vaporizer(R580), Inlet/Outlet End Caps(R580-IO) and tube accessory kit are included in the R540 anesthesia machine. Different air supply, induction chambers, small animal masks, filter canisters and gas evacuation apparatuses can be selected to form the experimental platform for anesthesia recovery of different animals

Suitable for rats, mice, hamsters, guinea pigs, rabbits, cats and other animals within 7kg

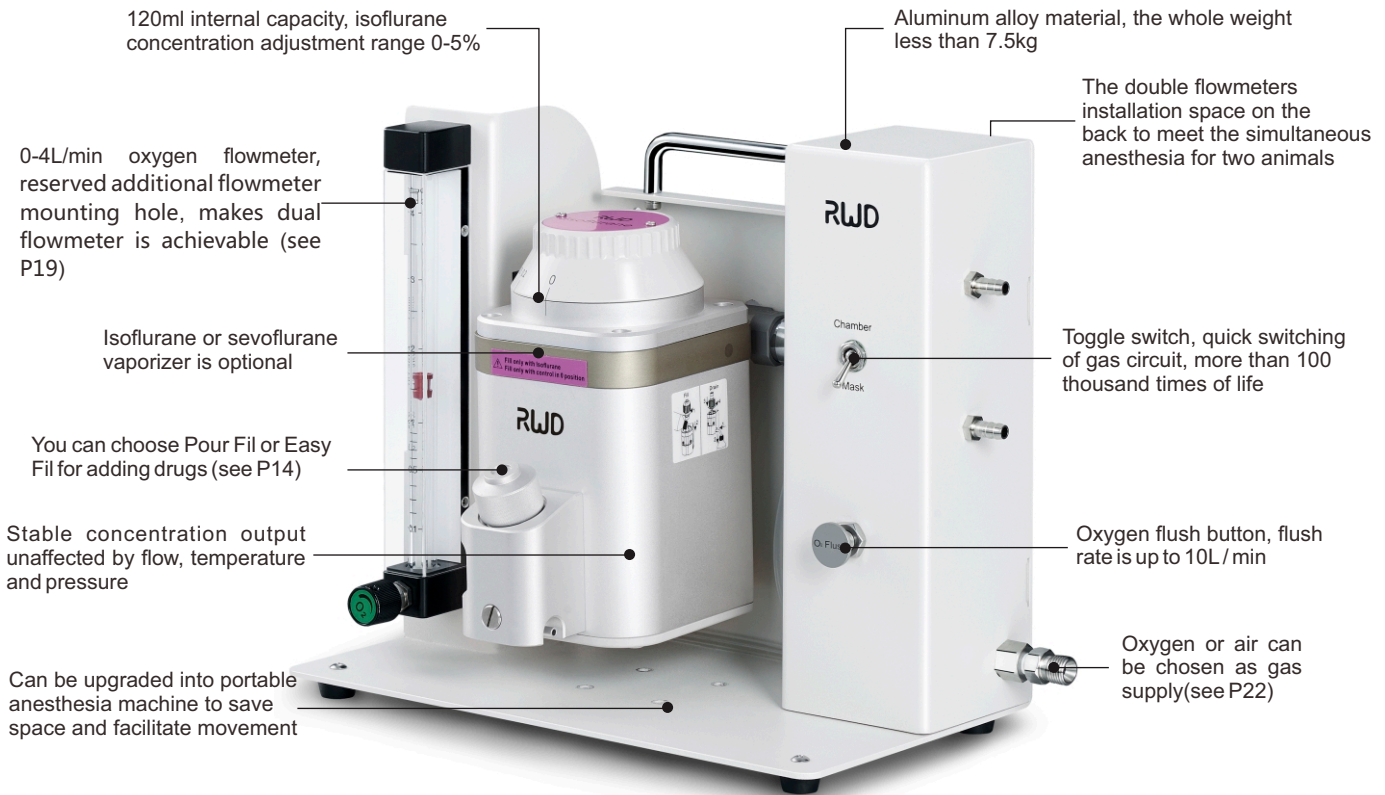

Note: Oxygen flush is commonly used to flush anesthetic gas to the filter canister before opening the induction chamber to avoid the escape of anesthetic gas into the personnel

Order Information:

Cat No.	Product Description	Configuration	Remarks
R540IE	Enhanced Small Animal Anesthesia Device-Isoflurane,Easy Fil	●	
R540IP	Enhanced Small Animal Anesthesia Device-Isoflurane,Pour Fil	●	Contains connection tubing for mask induction chamber
R540SE	Enhanced Small Animal Anesthesia Device-Sevoflurane,Easy Fil	●	
R540SP	Enhanced Small Animal Anesthesia Device-Sevoflurane,Pour Fil	●	
R510-S	Portable Stand for R510 Series Anesthesia Machine	○	Upgrade to portable type
R600-N2O-4E	Expansion Assembly with 0.1-4LPM N ₂ O Flowmeter	○	
R510-09	2ch Flowmeter for Anesthesia, range 0.1-1.0L/min each	○	Simultaneous anesthesia for two animals (adjusted independently)

Note: ●Recommended configuration ○ Optional