

RWD



Veterinary Oxygen Concentrator



One for Multi-use

supply oxygen for anesthesia, postoperative recovery, ICU, first aid and etc.



Cost-effective

Using only room air, it draws in and generates medical grade oxygen (USP93%), never run out of oxygen.



Quart-filtration System

Ceramic fiber • Intake filter cloth • Build-in filter • Humidifier bottle 99.999% bacteria, viruses and odor can be removed effectively.



Time Preset

Support presetting oxygen producing time, and will automatically stop working when time up.



Safe

Alarming Low Oxygen Concentration and Abnormal Situation, remind vet to handle on time.



Quiet


Quiet design, offers you a comfortable working environment.

One for Multi-use

ROC-5A/ROC-8A Veterinary Oxygen Concentrator is designed specifically for veterinary use. The high purity oxygen is separated from air by molecular sieve pressure swing adsorption (PSA) technology.

1.

Directly supply oxygen for animals

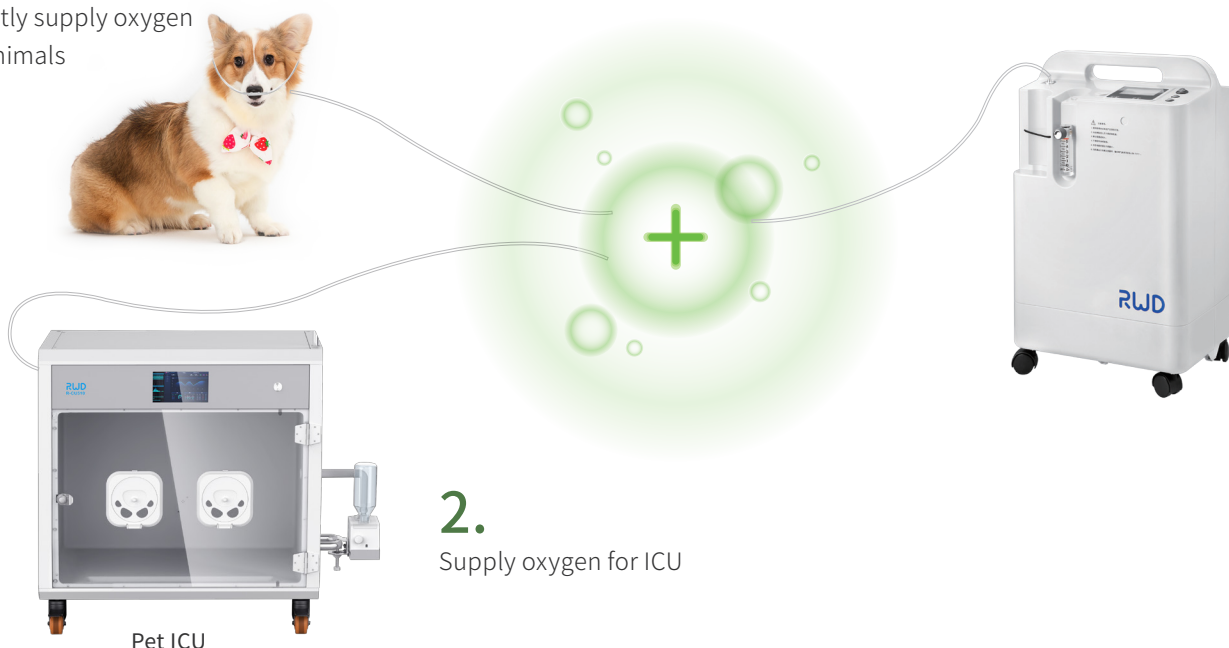

2.

Supply oxygen for ICU

Specification

Cat No.	ROC-5A	ROC-8A
Flow Rate	0-5L/min	0-8L/min
Net Weight	22kg	24kg
Power Rating	320VA	550VA
Veterinary Oxygen concentration	93%±3%	
Working voltage	~230 V/50Hz; ~110V/60Hz	
Outlet pressure	40-70kpa	
Dimension	300*360*600mm	
Noise	≤43dB	
Oxygen Tubing	2m*2	

3.

Supply oxygen for anesthesia



RWD Life Science Co.,Ltd

Add: 19-20/F, Building 7A, Shenzhen International Innovation Valley, Dashi 1st Road, Nanshan District, Shenzhen, Guangdong, P.R. China **E-mail:** rwd@rwdstco.com

RWD Life Science Inc.

Tel: (858)900-5879 **Support:** service@rwdls.com **Web:** www.rwdstco.com